
iRest

**iRest Supearl +
A336 MASSAGE CHAIR
OPERATION INSTRUCTIONS**

SAVE THESE INSTRUCTIONS

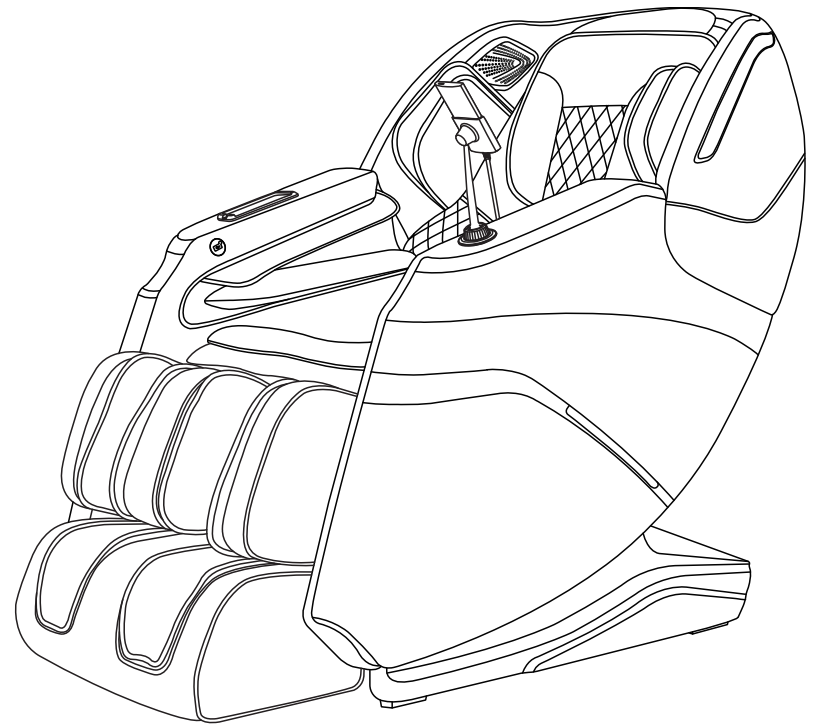Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our company reserves the right to change and explain the design of the product. If there is some difference between the product function in market circulation and the description in the manual, it is a normal update of the product, and our company will not inform you otherwise. The specific product function and color shall prevail in kind.

CONTENTS

Before Use

Safety Precautions	02-04
Installation Parts	05-06
Place Massage Chair	07-08
Before Use	09-10
Product Structure	11-14

Usage Method

Operation Instructions	15-42
------------------------	-------

Maintain

Cleaning and Maintenance	43
--------------------------	----

Important

Troubleshooting	44
Product Specifications	45

Safety Precautions must be observed

- Please read all relevant instructions carefully before use.
- The notes marked here will help you to use the product safely and correctly.
- In order to clearly indicate the extent of harm and damage, precautions are divided into **"Warning"** and **"Caution"** that may be caused by incorrect use. Both of them are important for safety and must be followed.

 Warning	Indicates that serious injury or death may result from improper use.
 Caution	Indicates that misuse may cause minor injury or property damage.

- Please pay attention to the symbol before each safety item.

 Forbidden	Absolute prohibition
--	----------------------

The following populations should not use this product

 Forbidden	■ The device is not intended for use by persons (including children) with physical disabilities, sensory or neurological impairments or who lack experience and common sense, unless they are supervised and instructed in the use of the device by those responsible for their safety. Children must be supervised to ensure they do not play with instruments. ■ Patients with osteoporosis. ■ Patients with heart disease or wearing electronic medical devices such as a pacemaker. ■ People who have a fever. ■ Pregnant women or women on their menstrual period. ■ Patients who are injured or have a skin disease. ■ When unsupervised, any children are under the age of 14 and those unconscious are not allowed to use this product. ■ People who are told by their physicians to have rest or who feel unwell. ■ People whose bodies are wet are forbidden to use this product. ■ This product has a hot surface. Persons insensitive to heat must pay attention when using this product. ■ If the supply cord is damaged, to avoid danger, it must be replaced by a professional from the manufacturer, its maintenance department or similar department
--	---

Safety Precautions must be observed

Important safety warnings

 Warnings	■ Do not allow children to touch any movable part of this product. ■ Use the well-grounded power supply unit suitable for this product. ■ Pull out the plug after use or before cleaning to avoid any physical injury or damage to the product. ■ Operate this product as instructed in this manual. ■ Do not use any parts or accessories which are not recommended. ■ Do not use this product outdoors. ■ Please read the Operation Instructions carefully before operation. ■ Never use this product for any other purpose not listed herein. ■ 20 minutes of usage each time is recommended. ■ Do not use this product if the cover or leather is broken or damaged. ■ Do not drop anything into this product. ■ Never fall asleep while using this product. ■ Do not use this product if you are drunk or feel unwell. ■ Do not use this product within one hour after a meal. ■ Do not use this product with excessive force to avoid any injury.
---	---

Safety precautions

 Caution	■ Check whether voltage is consistent with the specification of this product. ■ Never insert or pull out the plug with a wet hand. ■ Do not let water flow into this product to avoid electric shock or cause damage to this product. ■ Do not pull the power cord when plugging or unplugging it. Rough handling is forbidden. ■ Do not damage the wire or modify the circuit of this product. ■ Do not clean live parts of this product with a wet cloth such as a power switch and plug. ■ Be away from this product in case of power failure to avoid any injury if power is restored suddenly. ■ Stop using this product at once while it works abnormal, and consult your local supplier. ■ Stop using this product if you feel unwell, and consult your health care practitioner.
--	--

Storage conditions

 Caution	■ Storage room temperature: 5 °C to 40 °C. ■ Relative humidity (RH) between 20% and 80%; free of corrosive gases, and a well-ventilated room.
--	--

Product service and maintenance	
 Caution	■ This product should be maintained by the service center designated by the manufacturer only. No user should disassemble or maintain it without permission. ■ Do not forget to turn off the main power switch after usage. ■ Do not use this product if the power socket becomes loose. ■ If this product will be left unused for a long period of time, it is necessary to roll up the power cord of this product and store this product in a dry and dust-free environment. ■ Do not store this product at a high temperature or near open flame. Avoid prolonged exposure to direct sunshine. ■ Please clean this product with dry cloth. Never use thinner, benzene or alcohol. ■ The mechanical components of this product are specially designed and manufactured, so no special maintenance is needed. ■ Do not use a sharp object to stab at this product. ■ Do not roll or pull this product on any uneven ground. Instead, it should be lifted before moving. ■ Please use this product intermittently. Do not keep it working continuously for a long period of time.





Operating environment	
 Caution	■ Do not use this product in a very moist environment such as the bathroom. ■ Stop using this product immediately when the ambient temperature changes sharply. ■ Do not use this product in any heavily-dusty or caustic environment. ■ Do not use this product in a small space or a place without good ventilation.

Solutions to ordinary malfunctions	
 Caution	■ If the motor sounds during operation of this product, this is the normal operation sound. ■ If the controller cannot work normally, check whether the plug and the power socket are connected reliably, and whether the power switch is turned on. ■ If the rated working time is over, this product will automatically activate the power-off button; If this product keeps operating for a long period of time, the temperature protector will stop this product automatically. The product can be used again after half an hour.

Calf installation instructions

1. Insert the wire of the calf and the air pipe into the front socket of the steel frame assembly (Figure 1).
2. Then push the connected wire and air pipe into the hole in the backrest frame (Figure 2).
3. After the wire docking is completed, open the calf fixing covers on both sides of the steel frame assembly (Figure 3).
4. Lift the calf into the notch, close the fixing covers, and secure it in place with 2 pieces screws (Figure 4).





Note: this step can be ignored for all-in-one installation.

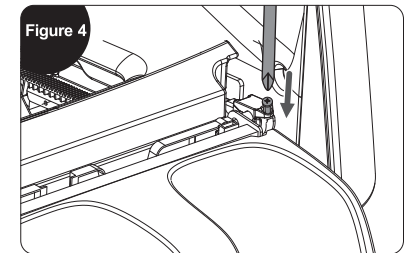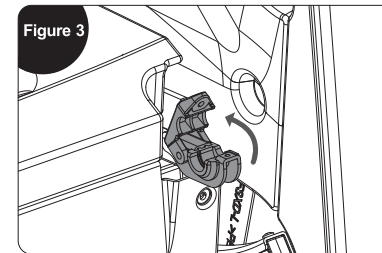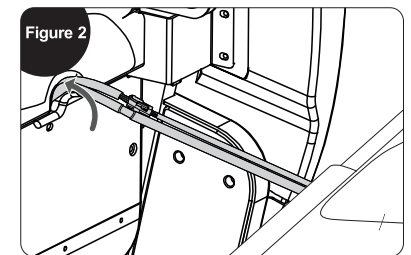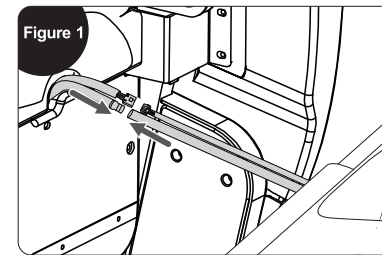Installation parts

Installation instructions for manual controller aerial plug

Non-fixed aerial plug installation method: insert the aerial plug (Fig. 1 Fig. 3) tightly after aligning it with the arrow mark of the aerial plug (Fig. 2 Fig. 4), and then lock the nut in a clockwise direction.

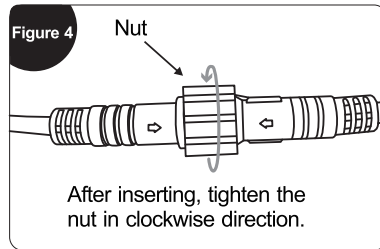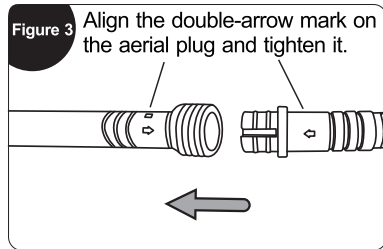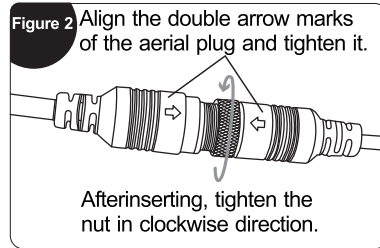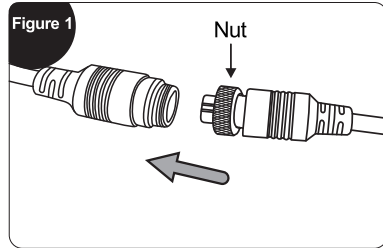
Place massage chair

Installation placement

Make sure there is enough space for the massage chair can be tilted at least 15 cm.

Do not expose the massage chair directly to high temperature environment due to sunlight (such as: in front of heating equipment), and place it to where can avoid high temperature and direct sunlight. It is recommended to lay a soft cushion for long-term use in a designated location. When laying the cushion, the size of the cushion should be enough to cover. Keep the range where the machine touches the floor and the range where your feet touch the floor.

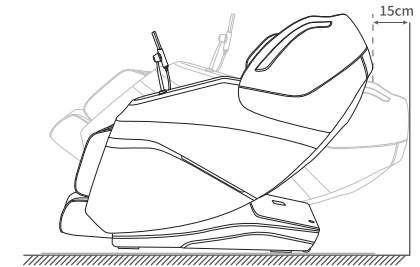
Product lift instruction

Require team work to lift and move the chair to avoid injury.

1. Please lift the massage chair to move it on vulnerable ground such as wooden floor.
2. When moving on a similar wooden floor, it is recommended to lift it. Avoid damage to the casters. When the floor is lifted and moved to the designated position, be careful to keep your hands out of the way to avoid injury, and be careful of your feet being crushed. Do not let go of the unit until the unit is completely flat, the floor may be damaged, therefore, it is recommended to place mats, etc. on the floor.

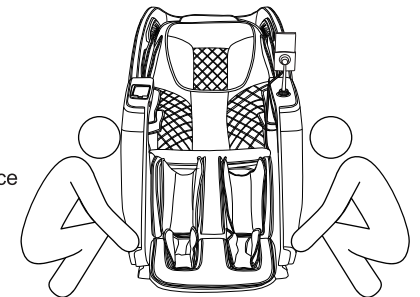
Place massage chair

Move method

Use the rear casters to move the massage chair. Move the chair by tilting it down and backward about 45° as shown in the diagram.



NOTE: Excessive force to tilt will cause the chair to tip over completely and product damage.

1. Before moving, please shut down the machine properly to reset the whole machine, and then cut off the power.
2. Please keep the chair vacant when moving. (People, pets or clutter will not be able to stay on the chair.)
3. When moving, do not cross a drop barrier higher than 2mm or a gap greater than 5mm.
4. When using casters to move this product, single continuous movement is not supported; if the distance exceeds 50 meters, please make this product move slowly and evenly.

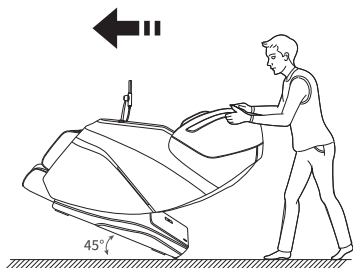How to disassemble and use headrest cushion and backrest cushion


Using a headrest cushion can reduce the intensity of kneading massage on the neck and shoulders, and you can decide whether to use a headrest cushion according to your needs (recommended to use a headrest cushion). The backrest cushion and the backrest are connected by a zipper (1), and the headrest cushion and the backrest cushion are connected by a zipper (2).

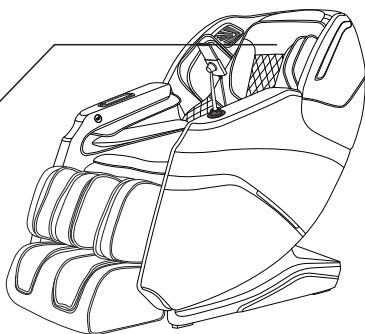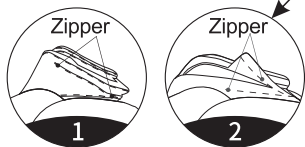Before use


Check around

Make sure that there are no people, pets and other objects around the unit.

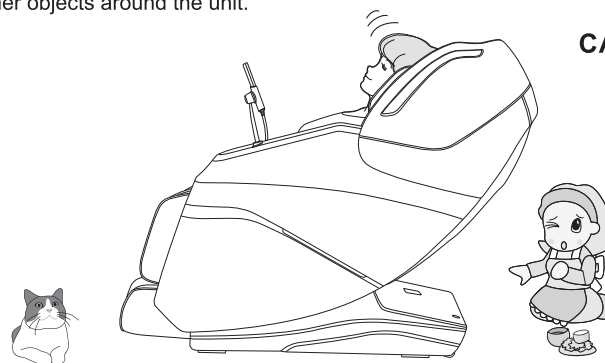Make sure there are no foreign objects in the gaps of the machine

Before massaging, make sure that no foreign objects are caught in the legs, feet, etc. Under no circumstances allow hands, feet, and head to penetrate into the area between the lower leg mechanism and main body, or the area between the armrest and the capsule. Under no circumstances should hands, feet, or head sink under the fabric or leather covering of the massage area.

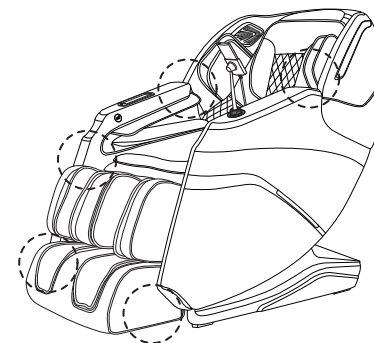
Before Use


Turn on the power

**Caution**


Plug in power cord

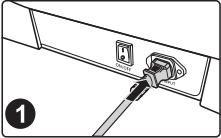Schematic diagram of power on

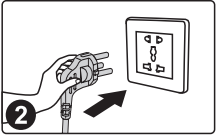Turn on the power switch ("I" position)

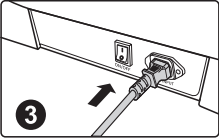**Warn**

■ Please check whether the power cord and power plug are damaged before use.

■ Please supply power to this product in accordance with the working power requirements of this product's calibration.

■ It is strictly forbidden to use a power supply other than the nominal power supply of this product.

■ It is strictly forbidden to connect various power conversion devices to supply power to this product without authorization.

10


Product Structure

Parts description


11

Product Structure

Function distribution

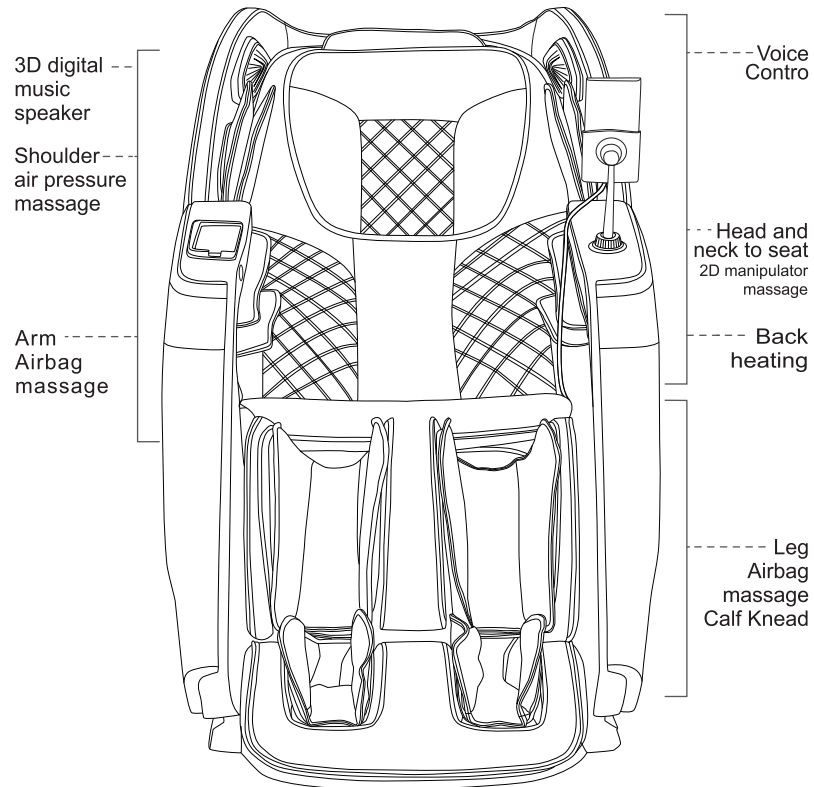
Product Structure

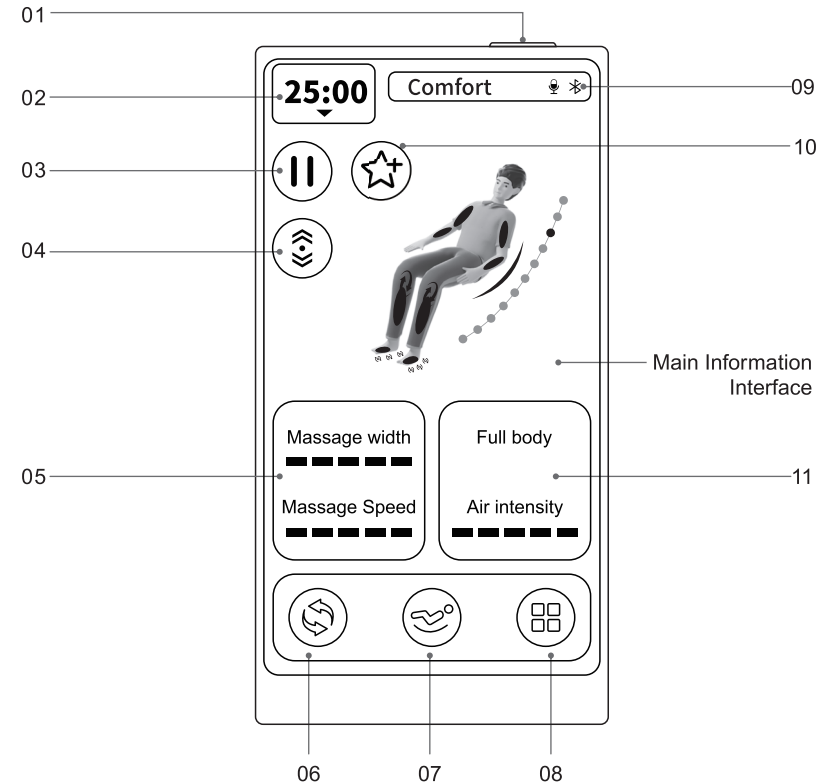01. On/Off key: Long press to turn on or off the massage chair.
02. Time adjustment key: Adjust the massage time.
03. Pause key: Pause all operations/ resume.
04. Shoulder position adjustment key: Enter into the shoulder position adjustment.
05. Massage Width/Speed Key: Enter the massage width/speed interface to adjust the massage width (3 adjustable levels) or massage speed (5 adjustable levels).
06. Auto key: Enter the automatic massage course selection.
07. Sitting posture adjustment key: Enter into the sitting position adjustment.
08. Menu key: Enter into the menu interface can choose automatic course, manuacourse, sitting angle, auxiliary function, memory, settings and other functions.
09. Status display bar: Display the current massage function status.
10. Save function key: Store the current massage program.
11. Air Pressure Massage/Intensity Key: Enter the air pressure massage selection interface to choose the corresponding air pressure massage and adjust the air pressure intensity (5 adjustable levels).

Product Structure

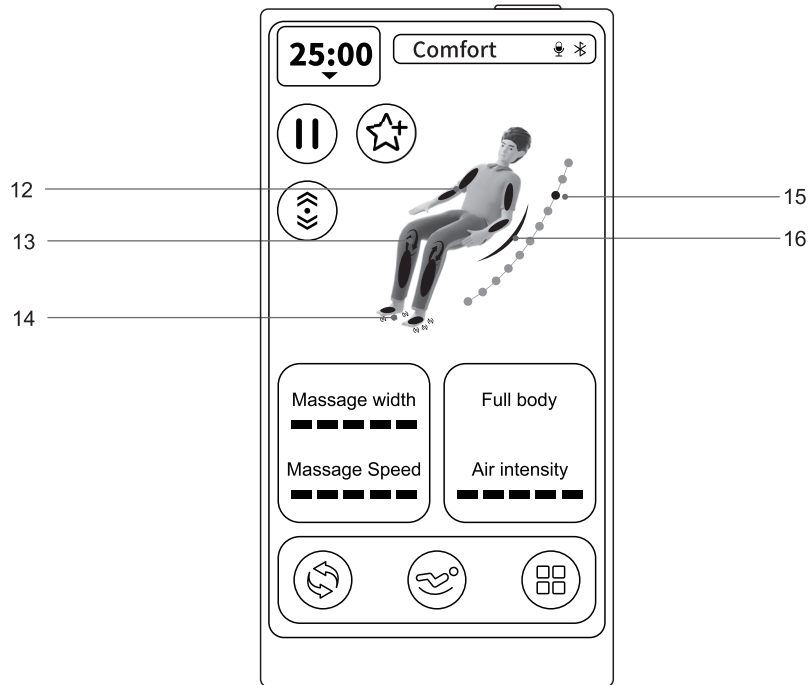
- 12、Air pressure display: Display the current air pressure massage position.
- 13、Calf knead: display the calf knead situation.
- 14、Foot roller: display the Foot roller situation.
- 15、Massage hand display: Display the current hand massage position.
- 16、Back heating display: Display the current back heating status.

Operation Instructions


I 、Start massage

1. Press and hold the power key at the top of the manual controller to turn on the massage chair. Then, wait for the manual controller system to boot up, and the system will automatically enter the "Automatic Recline" selection interface.

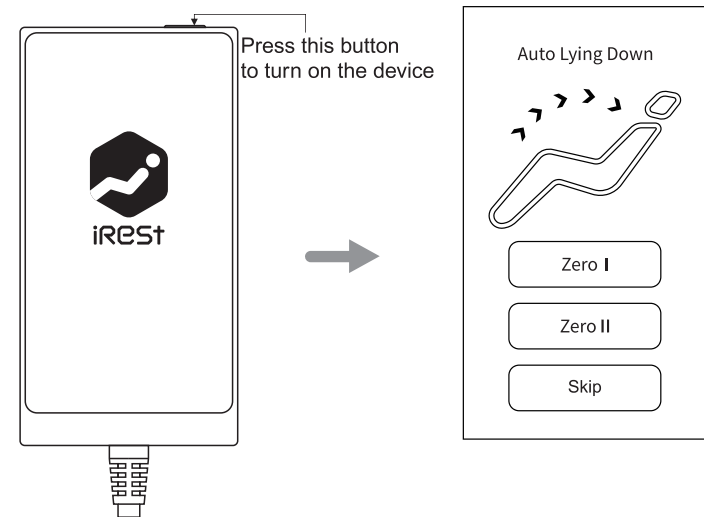
Operation Instructions

- (1).The angle of the massage chair is automatically adjusted after entering the automatic lying interface, and you can also manually select other zero-gravity angle, "Zero-gravity I or Zero-gravity II". When the angle adjustment is completed, the body shape detection will start automatically. If you do not need to adjust, you can also directly press the "skip" button, the massage chair will keep the current sitting position and proceed to the next step of body shape detection.
- (2).Enter the body shape detection interface, the massage chair will automatically detect the user's body shape. After the detection is completed, the massage ball will be positioned to the shoulder automatically, and enter the shoulder position adjustment interface. If you do not need to detect, you can also press "skip" key to use the system's default body shape and start the massage directly.

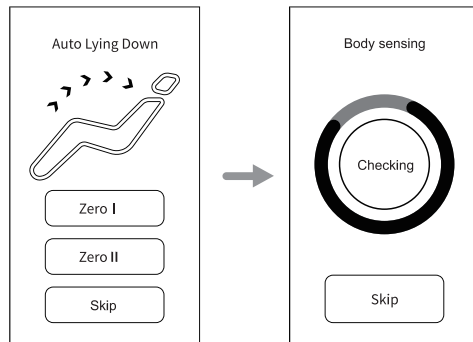
- (3).After entering the shoulder position adjustment interface, the manual controller will prompt you to manually fine-tune the shoulder position. As shown in the figure below, if the massage rollers stop at your shoulder position, no adjustment is needed—simply click "Confirm" to start the massage. If the massage rollers are not at your shoulder position, click the up and down arrows to adjust the rollers to the shoulder position (11 adjustable levels available). Once adjusted, click "Confirm" to start the massage. Note: If no action is taken for more than 10 seconds, the system will default the current roller position as your shoulder position and start the massage directly.

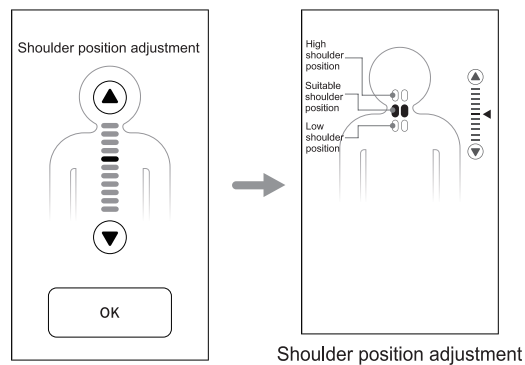Operation Instructions


- 2.Pause: On the home page, as shown in the figure below, click the "Pause" icon during program operation, and the massage chair will enter a paused state, stopping all massage functions. Click "Resume Massage" to restore the previous massage state. Note: In the paused state, except for powering off, other massage functions cannot be used.

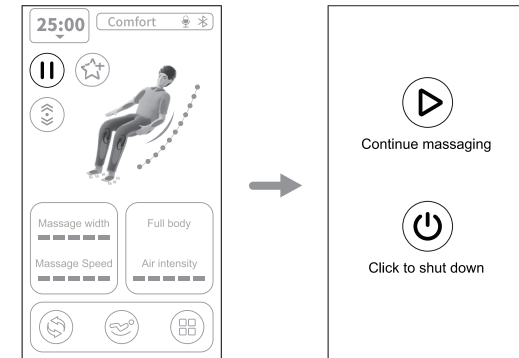Home Page-Pause/Running Diagram

- 3.Memory function: In the home page, as shown in the figure below, choose to press the "Save" icon to enter the save function selection interface, and then choose to store or recall according to the needs of massage. Choose to press the "Save" button, system will record the current state of the massage chair to run the program, choose to press the "Use" button, system will run the previously saved program. There are three sets of save modes for "Save" and "Use". Press the "Back" button to return to the home page.


Note: This function can also be operated from the main menu interface.

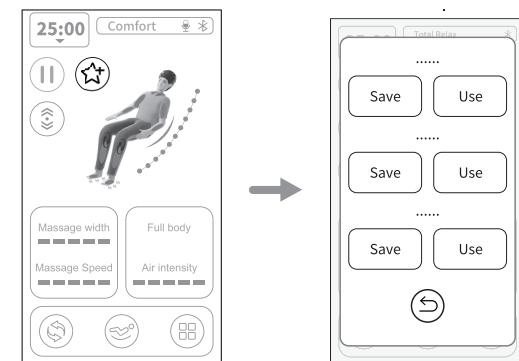Home Page -Save Diagram


Operation Instructions

4.Shoulder position adjustment: in the home page, as shown in the figure below, press the“**shoulder position adjustment**” icon to enter the shoulder position adjustmentinterface, and then adjust the position of the massage hand according to the massageneeds. If the massage ball is not in the position of your shoulder, press the up and downarrows to adjust the massage ball to the position of the shoulder (a total of 1-11 levelscan be adjusted), then press to confirm the selection and return to the home page.

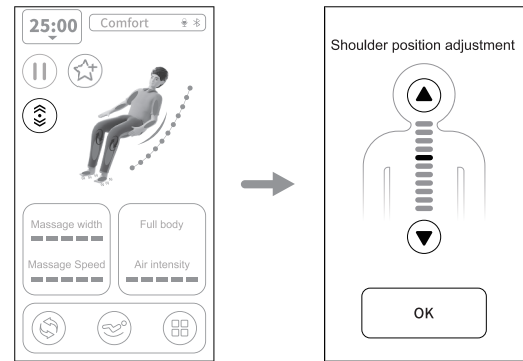Home Page-Shoulder Position Adjustment Diagram

5.Status Display: On the home page, as shown in the figure below, the current message function status information is displayed, such as Bluetooth, and other status details.

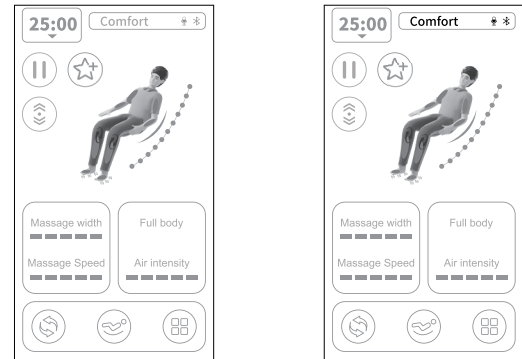
Home Page
Program Running Diagram


Home Page
Technique Display Diagram

Operation Instructions

6.Message Width Adjustment: On the home page, as shown in the figure below, select and click the "Message Width/Speed" key to enter the message width adjustment interface. Then, according to the message needs, slide the "Message Width" slider left or right to adjust the message width. There are 5 adjustable levels for width, with level 1 indicating the narrowest width and level 5 indicating the widest width. Click the "Return" key to go back to the home page.

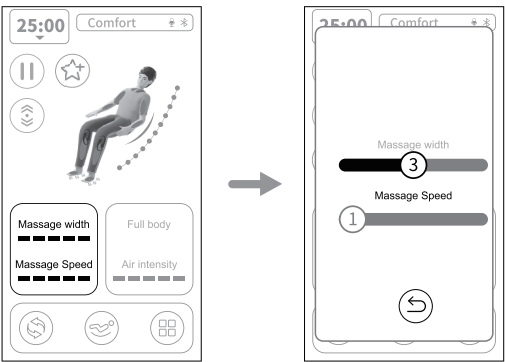
Home Page – Message Width Diagram

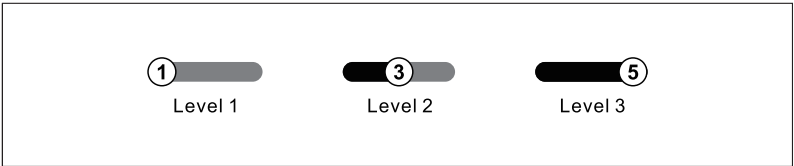
Button	Description
	Hold the slider key and slide to the right to increase the message width, adjustable from levels 1-3.
	Hold the slider key and slide to the left to decrease the message width, adjustable from levels 1-3.

Operation Instructions

7. Massage Speed: On the home page, as shown in the figure below, click the "Massage Width/Speed" key to enter the massage speed adjustment interface. Then, according to the massage needs, slide the "Massage Speed" slider left or right to adjust the massage speed, allowing free adjustment of the upper and lower massage arms' massage speed. There are 5 adjustable levels for speed, with level 1 indicating the slowest speed and level 5 indicating the fastest speed. Click the "Return" key to go back to the home page.
(speed is ok adjusted only on Manual program, Auto program cant adjust speed)

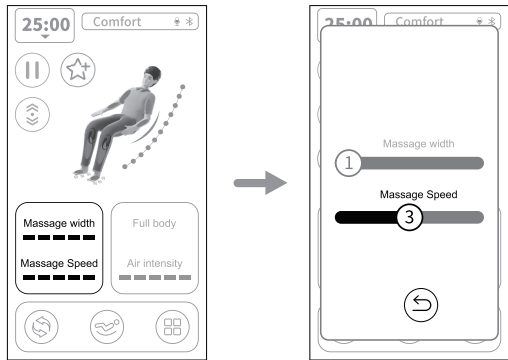
Home Screen - Massage Speed Diagram

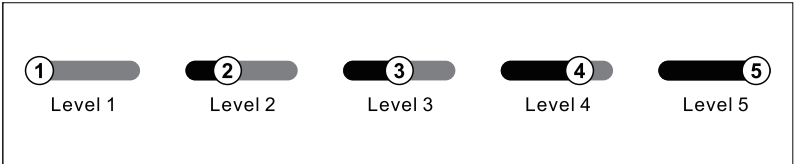Button	Description
	Hold the slider key and slide to the right to increase the massage speed, adjustable from levels 1-5
	Hold the slider key and slide to the left to decrease the massage speed, adjustable from levels 1-5

Operation Instructions

8. Air Pressure Mode: In the home page, as shown in the figure below, select the "AirPressure Massage /Air Pressure Intensity" key, enter the air pressure massageinterface, and then according to the massage needs to choose the corresponding ailpressure massage function. Air pressure modes including upper body, lower body,whole body, off and air pressure intensity adiustment, press off key is to turn off allthe air pressure function. Press the "Back" key to return to the home page.

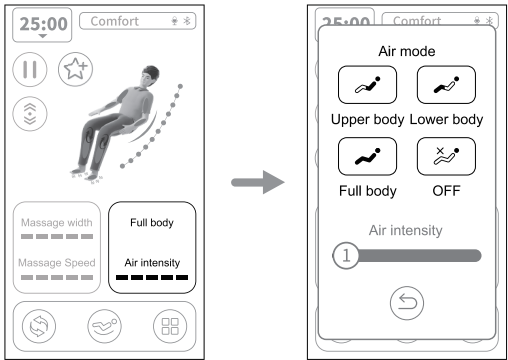
Home Page -Air Pressure Massage Diagram

Button	Description
	Press this key to turn off the air pressure function.
	Press this key to turn on the upper body air pressure massage function.
	Press this key to turn on the lower body air pressure massage function.
	Press this key to turn on the whole body air pressure massage function.


Air pressure display		
Upper body air pressure	Lower body air pressure	Whole body air pressure

Operation Instructions

(01).Air Pressure Intensity: in the home page , click the icon " air massage/air intensity" to enter to Air pressure massage operation surface. air massage mode including : Upper body ,Lower body , Full boady, OFF . Air intensity totally 5 levels , Slide the icon to adjust the intensity .

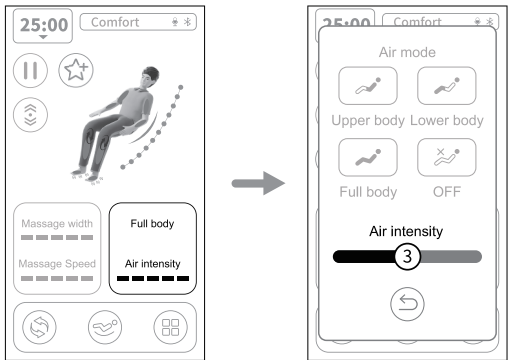Home Page -Air Pressure Intensity Diagram

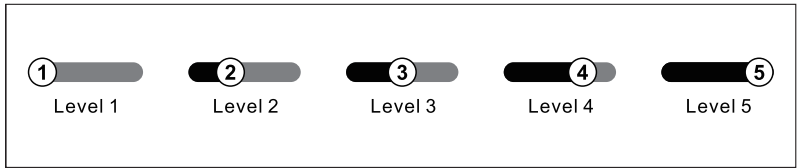
Button	Description
	Press and hold the slider to the right to increase the air pressure intensity, 1-5 levels of intensity can be adjusted.
	Press and hold the slider to move to the left to decrease the air pressure intensity, 1-5 levels of intensity can be adjusted.

Operation Instructions

II .Automatic Mode (please use this function after powering on the message chair)

Press the “Auto key” icon in the menu bar at the bottom of the screen as shown in the figure below ,the system will enter the automatic mode selection interface. Then according to the message needs to choose different automatic message course for massage. Note: This function can also be operated in the menu.

1. In the automatic mode interface, slide left or right to select from 12 automatic modes including Relaxing Massage, Comfortable Massage, Stretching massage, Back & Waist Massage, Neck & Shoulder Massage, Full Body Massage, Spine Care Massage, Fatigue Recovery, Meridian Treatment, Sitting Comfort, Soothing Sleep, and Leisure Health. Click the "Return" key to go back to the home page.

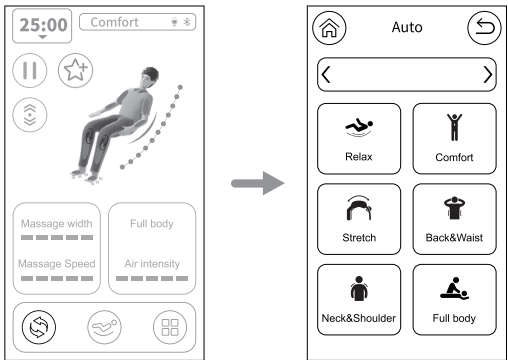Home Page -Automatic Function Selection Diagram


Operation Instructions

Overall Relaxation	
Function	Description
Comfort	Most comfortable massage experience.
Relax	Gentle and soothing deep massage to relieve body fatigue and promote relaxation.
Stretch	Frame and backrest frame rise and fall, and air pressure massage, and stretch or traction on the leg, remove fatigue, restore vitality.
Back & Waist	Focuses on the back and waist area.
Neck & Shoulder	Targets shoulders and neck for relief.
Full Body	Invigorating massage for the whole body.
Spine Care	Emphasizes spinal relaxation.
Fatigue Recovery	Relieves fatigue and improves microcirculation.
Meridian Relaxation	Promotes meridian balance and enhances the immune system.
Sitting Comfort	Reduces excess fat accumulation in hips and thighs.
Soothing Sleep	Improves sleep quality through acupoint massage.
Leisure Health	Combats sedentary lifestyle by improving body functions.

Operation Instructions

III. Sitting posture adjustment (please use this function after powering on the massage chair)

Press the **"Sitting posture"** icon on the menu bar at the bottom of the screen as shown in the figure below, enter the sitting posture adjustment interface, and then select the corresponding sitting posture adjustment function according to the massage needs, including linkage up, linkage down, backrest up, backrest down, calf up, calf down, zero gravity I, zero gravity II, reset and other sitting posture adjustment functions. Press the "Back" button to return to the home page. Note: This function can also be operated from the menu interface.

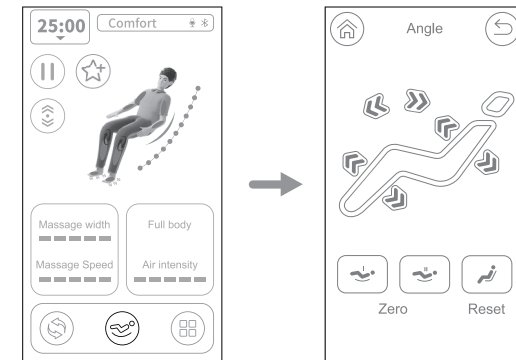Home Page -Sitting Posture Adjustment Diagram

1. Zero gravity adjustment

Button	Description
	Zero gravity I key: Automatically adjust the position of leg mechanism and backrest mechanism to zero gravity I position.
	Zero gravity II key: Automatically adjust the position of leg mechanism and backrest mechanism to zero gravity II position.

Operation Instructions

2.Calf angle adjustment

Button	Description
	Calf mechanism rising key: Press and hold this key, the calf mechanism rises slowly, release to stop rising.
	Calf mechanism descending key: Press and hold this key, the calf mechanism descends slowly, release to stop descending.

3.Backrest angle adjustment

Button	Description
	Backrest mechanism rising key: Press and hold this key, the backrest mechanism rises slowly, release to stop rising; when adjusting to the limit position, it will stop adjusting and emit a "be-be" sound.
	Backrest mechanism descending key: Press and hold this key, the backrest mechanism descends slowly, release to stop descending; when adjusting to the limit position, it will stop adjusting and emit a "be-be" sound.

4.Linkage angle adjustment

Button	Description
	Linkage rising key: Press and hold this key, the backrest mechanism rises slowly and the calf mechanism descends slowly at the same time, release to stop the backrest mechanism rising and the calf mechanism descending. when adjusting to the limit position, it will stop adjusting and emit a "di-di" sound.
	Linkage descending key: Press and hold this key, the backrest mechanism descends slowly and the calf mechanism rises slowly at the same time, release to stop the backrest mechanism descending and the calf mechanism rising; when adjusting to the limit position, it will stop adjusting and emit a "di-di" sound.


Operation Instructions

5.Reset

Button	Description
	Restore key: The angle of the massage chair is reset to its original position.

IV. Menu function (please use the function after powering on the massage chair)


Click the "Menu" icon in the menu bar at the bottom of the display screen, as shown in the figure below, to enter the menu function interface. Then, select the corresponding massage function according to your needs, including automatic functions, manual functions, posture adjustment, auxiliary functions, memory recall, settings, and other massage functions. Click the "Return" button to go back to the home page.

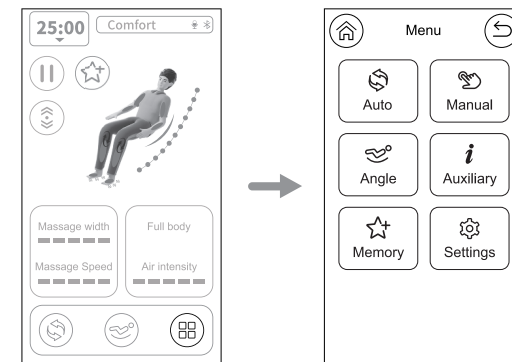Home Page -Menu Selection Diagram


1. Automatic Functions: After entering the menu interface, as shown in the figure below, click the "Automatic" icon to enter the automatic function interface. Then, choose the corresponding automatic function icon for massage according to your needs. Note: For detailed descriptions of automatic functions, refer to pages 22-23. Click the "Return" key to go back to the home page.

Operation Instructions

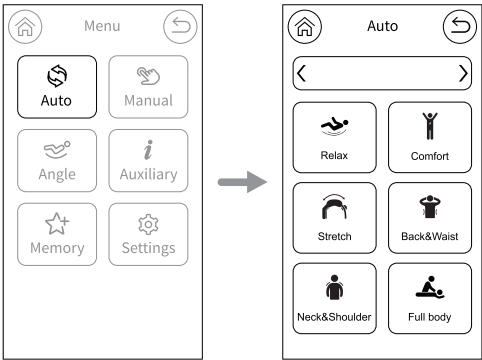Menu Interface -Automatic Function Selection Diagram


2.Manual function: After entering the menu interface, press the “Manual” icon as shown in the figure below, and then select the manual massage function according to the massage needs. Choose the “Back” key to return to the home page.

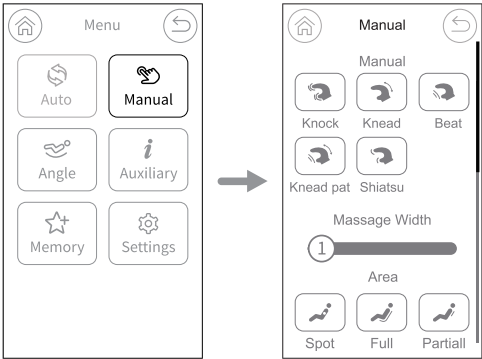Menu Interface-Manual Mode Diagram

01).Massge Techniques: After entering the manual interface, as shown in the figure below, click the "Manual" icon, then choose the corresponding massage technique icon according to the massage needs. Available techniques include, kneading, Beating, Kead pat , knocking and Shiatsu (5 techniques available)

Operation Instructions

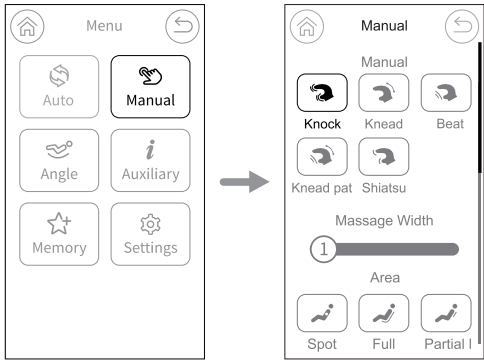
Menu Interface-Manual Mode Techniques Selection Diagram


Function	Massage roller adjustment	
Knock		Speed
Knead		Speed
Beat		Speed and width
Knead pat		Speed
Shiatsu		Width
		Adjustable speed in five levels Adjustable width in three levels

Operation Instructions

02).Massage Width: After entering the manual interface, as shown in the figure below, select and click the "Massage Width" slider button, then adjust the massage width by sliding left or right according to the massage needs. There are 3 adjustable levels for width, with level 1 indicating the narrowest width and level 3 indicating the widest width. Width adjustment is available for the following massage techniques: Lifting and Grasping, Kneading, Grasping, kneading, Beating, Kead pat , knocking and Shiatsu. Users must select a massage technique or massage area to enter manual mode before adjusting the width. Click the "Return" button to go back to the home screen.

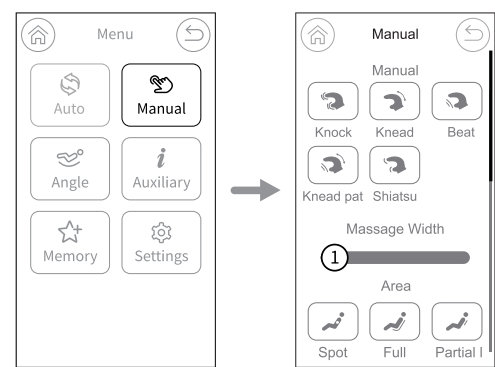
MenuInterface -Manual Mode Width Selection Diagram

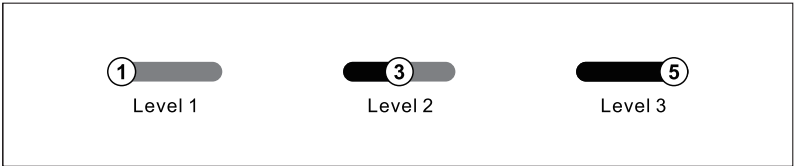Button	Description
	Press and hold the sliding key to slide to the right to widen the massagewidth, 1-3 levels of width can be adjusted.
	Press and hold the sliding key to slide to the left to narrow the massagewidth, 1-3 levels of width can be adjusted.

Operation Instructions

03).massage Area: in the manual surface , as show in the figure below choose the "Area"icon. there are totally 4 kinds of Massage Area for choose. Spot, Full, Partial I, Partial II. click the icon to start. press the "Back" button to turn to the home surface.

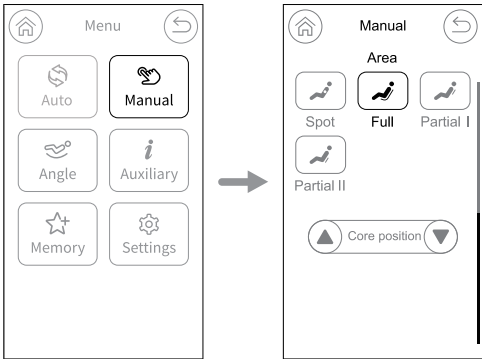
Menu Interface -Manual Mode Massage Area Selection Diagram

Button	Description
Spot	Fixed-point massage at fixed-point position press" Core position" to change the massage spot.
Full	Full back massage back and forth.
Partial I	Upper back massage back and forth.
Partial II	Lower back massage back and forth.

Operation Instructions

04).Message hand adjustment: After entering the manual interface, as shown in the figure below, select the core position according to the massage needs, press the "up arrow"or "down arrow" to fine-tune the massage hand up and down. Adjustment of the massage needs to be made in the manual mode of the fixed point press the "Back" key to return to the home page.

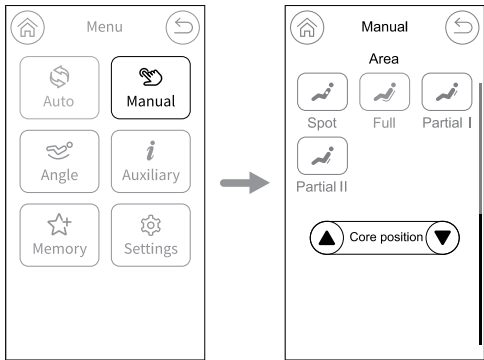Menu Interface -Manual Mode Message hand adjustment diagram


Message Hand Display Diagram

A 3D illustration of a person sitting in a chair. A series of black dots connected by a line represents the path of a massage hand along the person's back, starting from the lower back and moving upwards.

Button	Description
	Press and hold the upper key to move the massage hand upward, andrelease the key to stop moving the massage hand upward.
	Press and hold the lower key to move the massage hand downward, andrelease the key to stop the massage hand moving downward.


Operation Instructions

3.Sitting posture function:
Note: Seepages 27-29 for detailed description of the sitting posture functions. Press the "Backkey" to return to the home page.

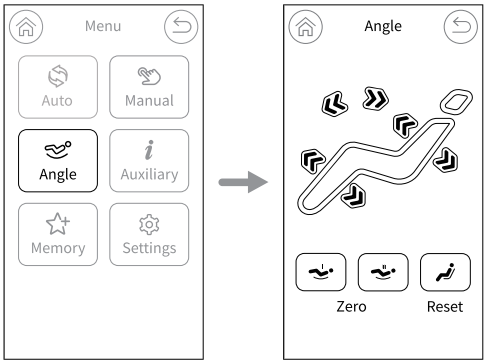Menu Interface -Sitting Posture Diagram


4.Auxiliary functions: After entering the menu interface, select and click the "auxiliary" icon button as shown in the following figure to enter the auxiliary function selection interface, and then select the corresponding auxiliary massage functions according to the massage needs, such as foot roller, back heating, calf kneadg and other massage functions. Click the "Back" button to return to the home page.

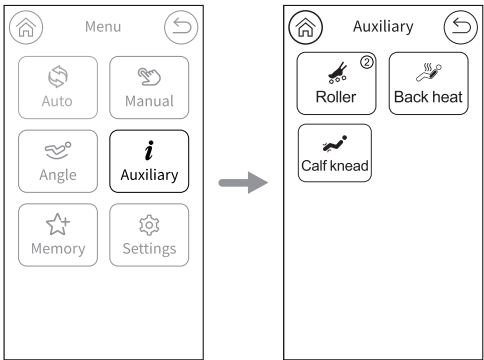Menu Interface -Auxiliary Function Diagram

Operation Instructions

(1).Foot roller: After entering the auxiliary interface, select and click the " roller" icon as shown in the following figure to turn on/off , and choose the intensity . there are 1-3 levels for choosing. Four state Loop switches.select and click "back" button to return to the home page interface.

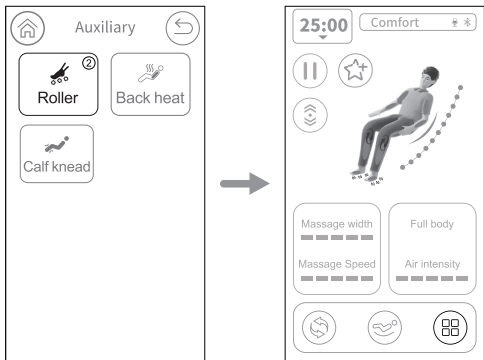
Foot roller diagram


Operation Instructions

(2).Back: After entering the auxiliary interface, as shown in the figure below, click the "Back Heating" icon, then turn on or off the back heating function according to the massage needs. The heating function will provide warmth after 3 minutes. Click the "Return" key to go back to the home page.


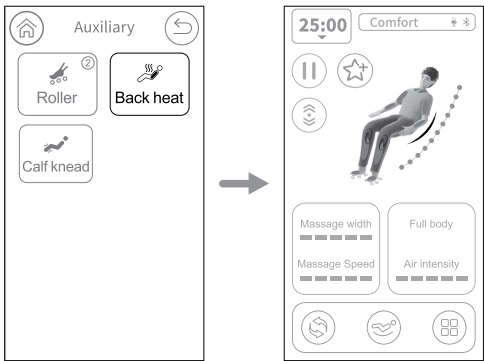
Home Heat display

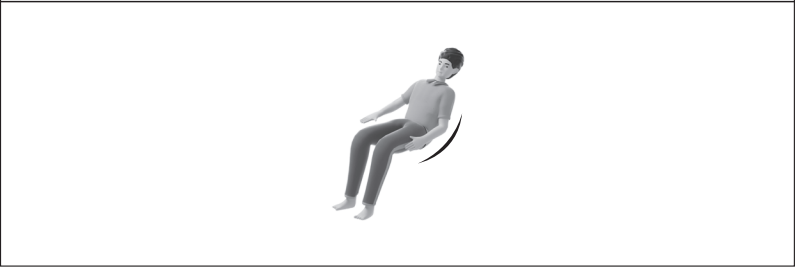
Operation Instructions

- (3).Calf knead: After entering the auxiliary interface, as shown in the figure below, click the "Calf knead" icon, then turn on or off the calf knead function according to the massage needs. Click the "Calf knead" key to go back to the home page.

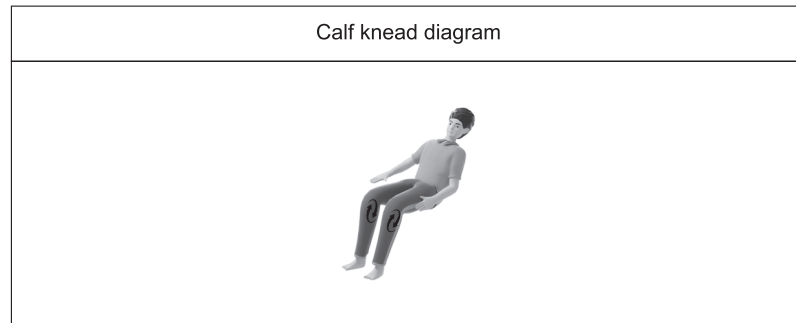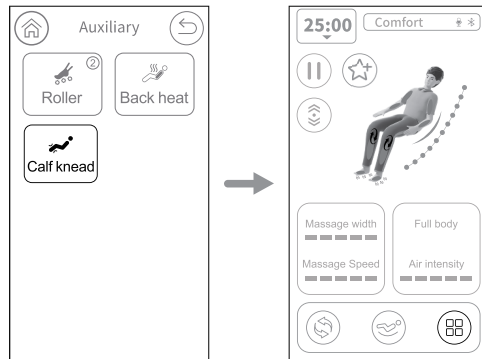
Operation Instructions

- 5.Memory Function: it is customized program by user self. easy for next time massage in the same status. Click "Memory" Icon, it will go to next surface. there are three Group "Save and Use". totally can store three customized programe. Save: click "Save" to stroge the program. it will save Current status of massage in use: autot or manual, chair lying angle, heating or not, intensity, ect Use: Click: "Use" to use the Saved massage.

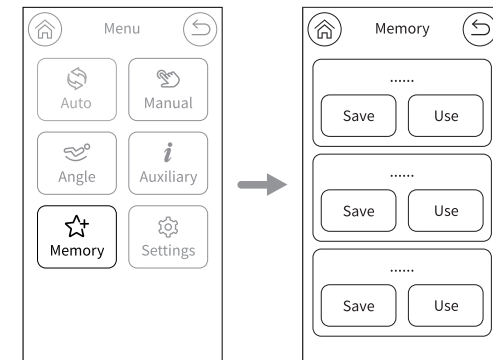Menu Interface -Storage Diagram


6. Setting function: in the Menu surface, there are "settings" icon. click this icon to find: Product information, Voice control, Bluetooth, Auto upright, Language. click "back" to back to home page.

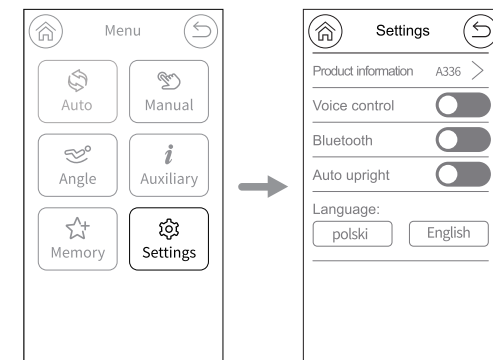Menu Interface -Setting Diagram

Operation Instructions

(1).Prdouct information: View Product information like product model, ID code, Software version.click "back" to back to home page.

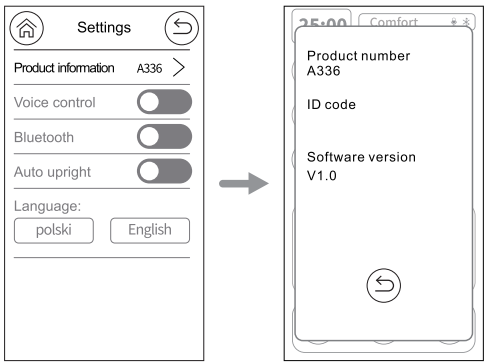
(2). Voice control: click to turn on/off voice control function. Click "back" to back to home page.

- 1)、 When the massage chair is turned on or in standby mode the sound "Hi Alice"or "Hey Alice"will be produced to enable the voice control function.
- 2)、 After such function is enabled successfully, the massage chair will reply "I am here" At this moment, once a corresponding voice control entry to be implemented is given within six seconds, the corresponding voice control entry to be implemented is given within six seconds, the corresponding voice control function can be performed after a reply is received from the massage chair if no correct voice control entry is given within six seconds, the voice control function automatically will be off automatically and needs to be enabled again.
- 3)、 Note: When using the voice control function, the speech speed should not be too fast, and the voice commands should be issued according to the voice control entries. After each successful use of the voice-activated entry, it is recommended to give out the voice command after 1-2 seconds.

Operation Instructions

Voice Control Entry		
Entry: "Hi Alice" or "Hey Alice". Answer: I am here.		
No.	Voice control entry	Answer entry
1	Massage on	Answer: Ok, massage on
2	Massage close	Answer: Ok, massage close
3	Comfort Massage	Answer: Ok, comfort massage
4	Relax massage	Answer: Ok, relax massage
5	Full body massage	Answer: Ok, full body massage
6	Neck and shoulder massage	Answer: Ok, neck and shoulder massage
7	Back and waist massage	Answer: Ok, back and waist massage
8	Stretch massage	Answer: Ok, stretch massage
9	Open the air pressure	Answer: Ok, open the air pressure
10	Close the air pressure	Answer: Ok, close the air pressure
11	Up the seat position	Answer: Ok, up the seat position
12	Down the seat position	Answer: Ok, down the seat position
13	Change the other mode	Answer: Ok, change the other mode
14	Go little down	Answer: Ok, go little down
15	Go little up	Answer: Ok, go little up
16	Fixed massage	Answer: Ok, Fixed massage
17	Stop adjustment	Answer: Ok, Stop adjustment

Operation Instructions

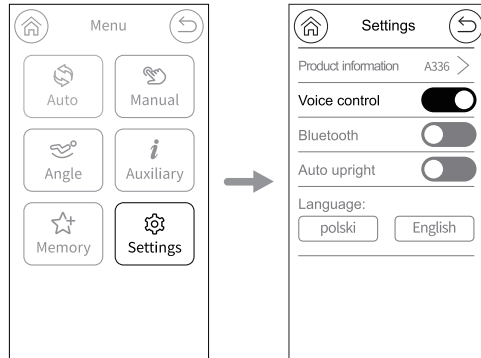
(3).Bluetooth: After entering the setting interface, press the **"Bluetooth"** key as shown in the figure below, and then turn on or off the Bluetooth function according to the massage needs. After turning on the Bluetooth function, you can pair and connect the audio source device with Bluetooth function (such as cell phone, MiD tablet PC, etc.) with the Bluetooth module of the massage chair, and then the music played by the audio source device can be transmitted wirelessly via Bluetooth to the sound system of the massage chair for music playback. Choose to press the "Back" key to return to the home page.

01) Enter the smart device main menu, Open "Settings" and select the "Bluetooth" options. Open "Bluetooth", and Bluetooth will automatically "Search for a New Bluetooth Device".


02) The smart device will search for the corresponding Bluetooth model (e.g. IMCM-XXXX), Click the "Match" button, if you are asked whether to match with the device, you just select "Match" to complete matching.

03) Open the music player of the smart device and select a song to play. Press the "Back" key to return to the home page.

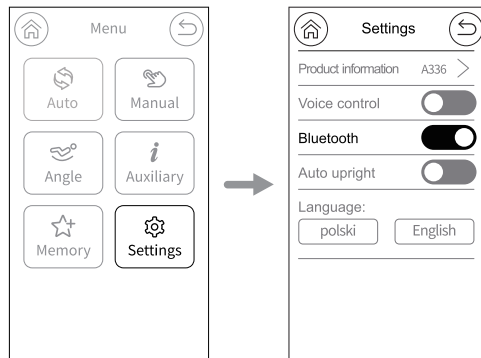
Operation Instructions

(4).Auto upright: it is a Default setting . Turn on the "Auto upright " the chair will Default automatically reset when massage finished. Turn off "Auto upright" ,the massage chair will default not reset after massage finished. please setting according to your Usage habit.

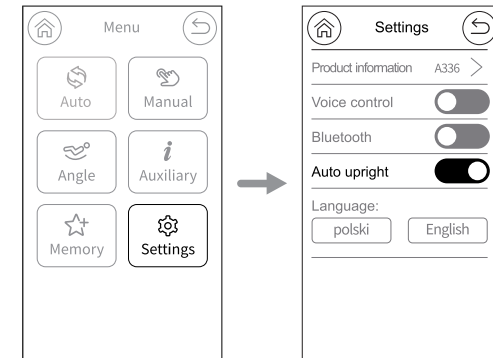
(5).Language Selection: After entering the settings interface, as shown in the figure below, set the system language according to the massage needs. The system offers 2 language options, with English as the default and Polish as an alternative. Click the "Return" key to go back to the home page.

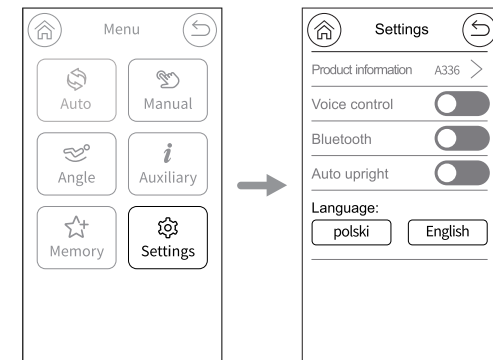
Operation Instructions

V、End massage in the massage process:


1. Long press the power button on the top of the hand control, all massage functions off, massage chair angle restore.
2. Press the pause key, then press the off key on the hand control interface, all massage functions off, massage chair angle restore.
3. Use the voice control command to power off, which will close all massage functions and reset the massage chair angle.
4. The massage time is over, all massage functions off, if the auto upright function is turned on, then the massage chair angle is restore, otherwise the massage chair angle remains unchanged.

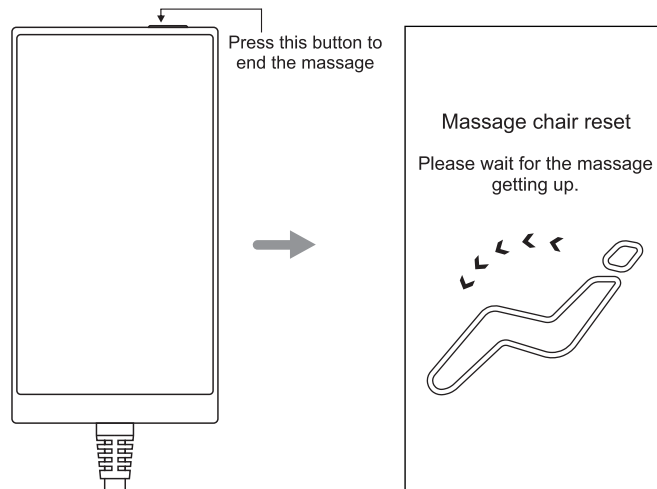
Turn off the power to the entire machine.

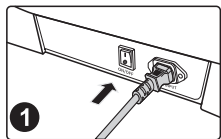Power off ("0" position)

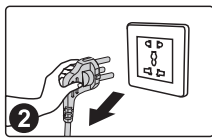Take off the power plug

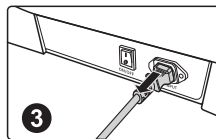Take off the product power cable.

Cleaning and Maintenance

Clean PU leather and plastic parts

Wipe the leather with a soft, dry cloth.

- Do not use chemicals such as thinner, gasoline, alcohol, etc. to clean the unit.
- When using commercially available leather care products (rags), please follow the instructions for use.
- If the PU leather is particularly dirty, wipe it by the following methods.

- ① Immerse a soft cloth in water or hot water containing 3% to 5% neutral detergent, then wring it out thoroughly.
- ② Use the above soft cloth to wipe the leather surface.
- ③ Wipe off the cleaning solution with a cloth soaked in clean water and thoroughly wrung out.
- ④ Wipe with a wrung soft cloth.
- ⑤ Let it air dry naturally.


- When it is difficult to remove stains, please immerse commercially available melamine foam in neutral detergent, and then wipe the machine with it.

- Do not use a hair dryer to dry the surface.

- As the color of the fabric will stain the surface of the PU leather, be careful when wiping the chair with denim and colored fabrics.

- Do not keep these parts in contact with plastic for a long time, as this will cause discoloration.

- PU leather may be discolored, so when using hair dye, please cover the touched part with a towel.

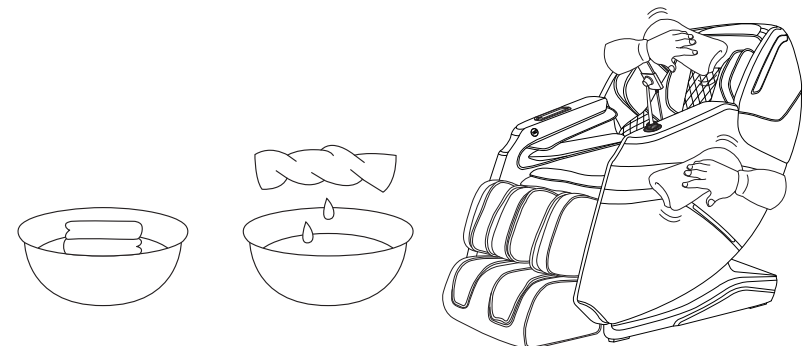Troubleshooting

Fault	When the calf mechanism or the backrest reaches a certain position, press the corresponding adjustment button and it will no longer change, and a continuous "beep, beep" sound will be emitted.
Solution	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.
Fault	The function does not work or some functions do not work after the machine is turned on.
Solution	After the machine is turned on by pressing the switch key, it is necessary to manually select an automatic massage program or other working state. If it is not operated, the machine will automatically shut down after 20 minutes. If the manual operation mode is selected, the functions that are not selected may be disabled, and the working mode needs to be manually set in the manual mode.
Fault	The machine makes a leather rubbing sound when adjusting its posture.
Solution	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.
Fault	The massage wheel did not reach the shoulders or neck.
Solution	If the head does not touch the head cushion or the back does not touch the backrest, the shoulder position may be lower than the actual position during the body shape detection process. Sit in the deepest position of the seat, rest your head on the headrest, and start over from the beginning.
Fault	The unit is damaged. The power cord or power plug is abnormally hot.
Solution	To prevent accidents, be sure to contact an authorized service center.

Please feel free to contact the designated service personnel for any questions or doubts about this product.

Product Specifications

Model: A336

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz ☐
220-240V~ 50Hz/60Hz ☐

Rated power: 140W

Rated time: 20min

Safety structure: Class I

Material description: PVC, PA, steel parts, electric and electron components

MATERIAL	PU	LEATHER	CLOTH	WOOD

SAVE THESE INSTRUCTIONS

Producer:

Zhejiang Haozhonghao Health Product Co., Ltd

No.18 XingLong Road, Furniture Garden, Wanquan Industrial Base, Pingyang, Zhejiang, China

e-mail: irest@irest.com

Importer and authorized representative in the EU:

iRest Europe Sp. z o.o.

Szewczeni 7/11

10-274 Olsztyn, Poland

e-mail: irest@irest-europe.com

Dear user, you can contribute to environmental protection. Be environmentally responsible and return your used equipment to the relevant collection point designated by your local administration. The materials used to manufacture the device may be valuable recyclable resources and, if disposed of incorrectly, can pose a risk to the natural environment. Remember - do not dispose of the device together with your household waste.

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